


LCP Mixing Station

Spend your time screening your protein, not mixing your LCP. Small, quiet and easy to use, the LCP Mixing Station provides consistent mixing results from sample to sample, user-to-user and day to day. Simply loading the syringe automatically sets the mixing volume so that only the sample volume is mixed.

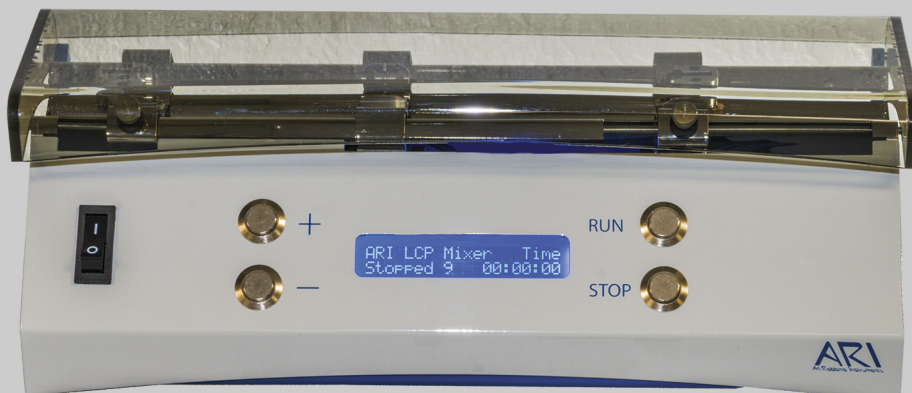

Features

- Eliminate user-to-user and day-to-day variation in LCP setups.
- Accommodates Hamilton style syringes and/or the ARI LCP dispensing syringe. Any volume syringe can be used provided both syringes are the same volume.
- Quick and easy to load.
- Loading the syringe automatically sets the mixing stroke so that the mix volume is equal to the sample volume in the syringe.
- Elapsed mixing time is displayed.
- On the fly mixing speed adjustments.
- Easy to clean polished stainless steel mixing chamber.

Specifications

Catalog #: 620-4000-00 Includes mixing coupling. Ferrules and syringes sold separately.
Size: 13.5 (34.3) x 5.5 (13.9) x 4.5 (11.4) inches (cm)
Weight: 5.5 lbs, 2.5 Kg
Electrical: 12 Vdc, 2 Amps
Environment: 4° to 70° C, Non-condensing



For more information, contact us at: Tel: (888)658-5300
www.artrobbins.com
info@artrobbins.com


LCP Handheld Dispenser

Handheld LCP Dispenser

Designed specifically for dispensing LCP at low volumes with 10 μ L syringe. It is also compatible with other syringe sizes. A perfect solution for a lab just getting started or a larger lab that just needs to do a few quick optimization setups.

The volume of the syringe can be divided into 50, 100 or 200 steps per dispense. With a 10 μ L, the dispense volume can be set to 50, 100, or 200 nL. This variable volume feature is a must have for crystal optimization.

Catalog: 620-4000-00 Handheld Dispenser
620-4101-00 10 μ L Syringe
620-4102-00 Syringe Needles (3x)

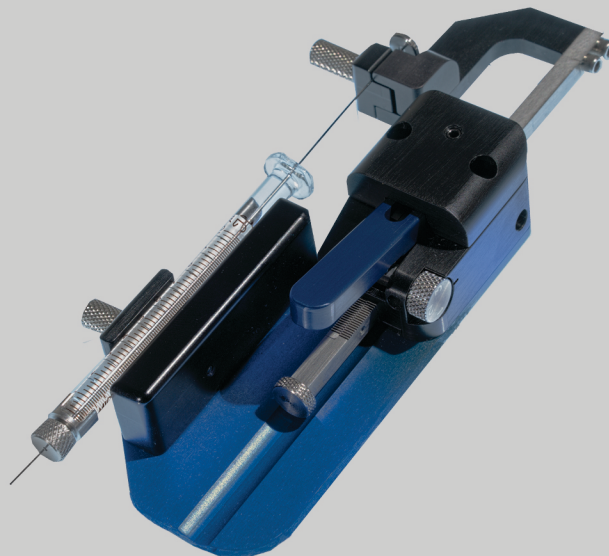


LCP Mixer/Dispenser Starter Kit

Get the LCP mixer and handheld dispenser as a kit and save money.

Catalog #: 620-4250-00

Includes:

- LCP Mixer (1x)
- LCP Handheld Dispenser (1x)
- 10 μ L Syringe (1x) + Needles (3x)
- 100 μ L Syringe (2x)
- ARI-ARI Ferrule (1x)
- ARI-Hamilton Ferrule (1x)

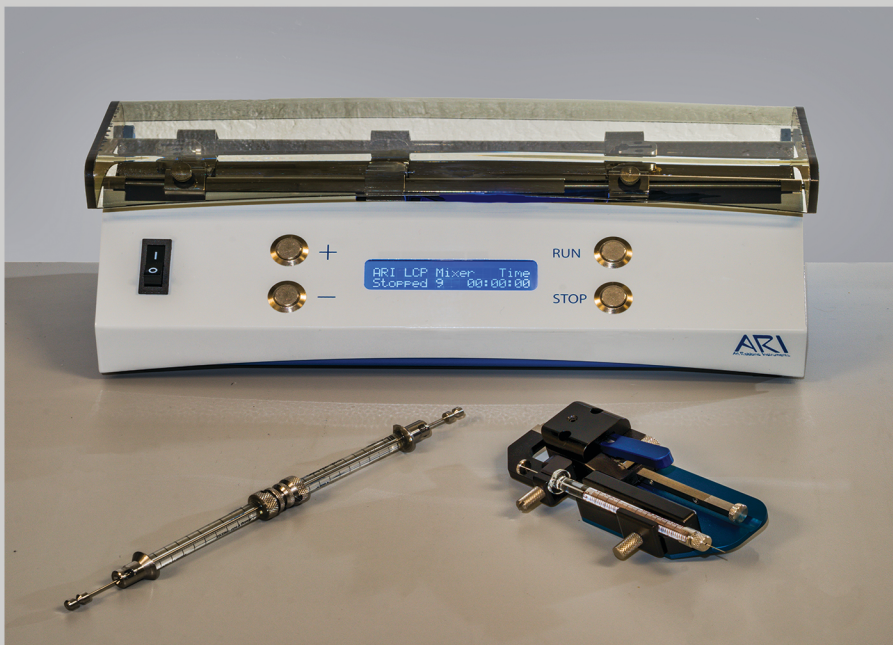

For more information, contact us at: Tel: (888)658-5300
www.artrobbins.com
info@artrobbins.com